



Healthward uk

Healthier your life

CD3q

FEATURES AND HIGHLIGHTS







Technical parameters

External size (LxW)	2030x1045 mm
■ Bed height	500-770 mm
 Mattress platform 	1705x900 mm
Back-rest adjustment	0-75°(±10°)
Knee-rest adjustment	0-35°(±10°)



Technical configuration

✓ Manual Cranks	3pcs
✓ 5" Hardshell Castors	4pcs
Aluminium Alloy Full Length Side Rail	1set
ABS Grace Bed Ends with Safe Lock	1set
✓ IV Pole	1pc
✓ IV Pole Prevision	4pcs
☑ Drainage Hook	2pcs





Design

Head foot board, with curve surrounded design, provide conforting feeling for patients.



Plastic

Environmental plastic can be sterilized in 100°C, tensile upto 30MPa Virgin ABS material.



Steels | Painting

Baosteel from Fortune Global 500 11 Process epoxy painting, ASTM testing anti-baterial, paint thickness 0.12 mm, brightness 60°, paint can resist 50 kg impact.



Welding | Chassis

Panasonic robotic ensure, 360° full smooth wedlding Solid bed base made of 5.0 mm stamped rotating plate.



Hand Control

Retractable hand crank, 100,000 times remaining intact Extreme protection device, sounds when tipped to the max position.



Platform

Stamping forming platform, smooth surface no welding scar, bearing 250KG load. Option resin X-ray translucent platform.



Guardrail

Sliding full side rails with distinctive design, elegant appearance; The lock bar is in conspicuous red,easy operation, sliding upper and down smoothly; Anti-pinch design, effectively protect the safety of children in use.



Head & Foot Board

Head and foot board can be reversable, extra height design.

ABS Grace Bed Ends with Safe Lock, automatic lock without any manual operation on lock, on/off removable mark on hook lock.



Castors

5" hardshell casters TPR tire no worn out after running 30KM, save anti-winding hard shell, united forming without bolts. Pass Dynamic test: Bearing 120kg run 30KM and pass obstacles 500 time.



Drainage Hooks

Mobile plastic drainage hooks.



IV pole

Stainless Steel Column with plastic hook.

















Design

Manufacturing

Maintenance



USA

Saikang Medical Technology, Inc. 17800 Castleton St Ste 180, City of Industry, 91748, California. Tel.: +1 949648 2776

United Kingdom

UK Healthward International Ltd. Dept. 111, 196 High Road Wood Green, N22 8HH, London. Tel.: +44 2032399738

Russia

Saikang Medical Russia LLC Truda 174, 454080 Chelyabinsk. Tel.: +7 (351) 214 5566 / 8811 / 1516

East Africa

Saikang Medical Solutions Ltd. Viraj Complex Unit 3, Mombasa Road, P.O. Box 14262-00800, Nairobi. Kenya. Tel.: +254 733 704 690

Middle East


Saikang Medical Supplies LLC X22-S03 England Cluster, International City, Dubai. United Arab Emirates. Tel.: +20 128 8962497

Healthward uk



No.1, 5th Floor, Building 5, No.360, Wuke West 5th Road •



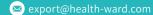

(www.health-ward.com







Chengdu Saikang Medical Equipment Co., Ltd. STOCK CODE: 870098

Company video

General catalogue