



Healthward uk

Healthier your life

Product data sheet



FEATURES AND HIGHLIGHTS

Technical parameters

External size (LxW):	(2040-2240)x1020 mm
Mattress platform:	(1820-2020)x820 mm
Safe working load:	250 kg
Height adjustment:	500-950 mm
Back-rest adjustment:	0-75° (±10°)
Knee-rest adjustment:	0-20° (±10°)
Reverse Trendelenburg:	0-18°
Trendelenburg:	0-18°
Angle indicator:	0-90°
 Voltage and frequency: 	100-240 V · 50/60 Hz
 Plug as required (hospital grade) 	

Technical configuration

☑ Linear LINAK Motor	6 pcs
☑ 6" Double Side Castors, Central Locking	4 pcs
☑ ABS Side Rail with Integrated Control Panel	1 set
☑ ABS Bed Ends One Step Removable	1 set
☑ Ø75mm Wall Bumper Wheel	4 pcs
☑ IV Pole	1 set
✓ IV Pole Prevision	4 pcs
☑ Drainage Hook	4 pcs
☑ 5th Directional Wheel Pedal	2 pcs
☑ Central Locking Pedal	2 pcs
X-ray Translucent Mattress Platform	1 set
☑ Back-up Battery	1 set
☑ Night Light	1 set
☑ Bed Push Handle for Emergency Transfer	1 set
✓ Manual CPR	
☑ Electric CPR	



























Optional

☐ Weighing System With LCD Screen
☐ Data Storage System
Mankey Pale

Design

Extra Long and height head foot board, with curve surrounded design, provide conforting feeling for patients.



Welding | Safety

Panasonic Robotic ensure, 360° full smooth wedlding. IEC606011 Safty Standard.



Electric motor

Linear motor, with UL/ROHS/EN standard. Noise under 50 dB, with lifetime 20,000 times movement.



Brake

Break spring keep qualified after 10,000times uninteruped break.



Bed surface

X-Ray Translucent Mattress Platform.



ABS bed ends

ABS bed ends One Step Removable, meet the first-aid requirements of high-end wards and ICU room.



Battery

Battery - 18times movement with 100kg bearing, 130times movement without bearing.



Bottom

The shelf at the bottom of the bed can hold patients items, like shoes and something others.



Steels | Painting

Baosteel from Fortune Global 500. 11 Process epoxy painting, ASTM testing anti-baterial, paint thickness 0.12mm, brightness 60°, paint can resist 50kg impact.



Plastic

Environmental plastic can be sterilized in 100°C, tensile upto 30MPa.



7 Inch touch screen

Replacing the orginal button screen, more convenient operation. There are two language options, chinese and english. Screen can store the patient and doctor's information, users manuals, weight record



Wheels

TPR tire no worn out after running 30 km, save anti-winding hard shell, united forming without bolts. Pass dynamic test: Bearing 120 kg, run 30 km and pass obstacles 500 time.



Side rail

Tuck away side rail with buffer Built-in controller and angle indicator.



Bumper wheels

Bumper wheel , ϕ 74mm, max. protect the bed and patient during transportation.



IV pole

Stainless Steel Column with plastic hook.



Lifting pole

Lifting pole - Bearing 765N force.

















Design

Manufacturing

Maintenance



USA

Saikang Medical Technology, Inc. 17800 Castleton St Ste 180, City of Industry, 91748, California. Tel.: +1 949648 2776

United Kingdom

UK Healthward International Ltd. Dept. 111, 196 High Road Wood Green, N22 8HH, London. Tel.: +44 2032399738

Russia

Saikang Medical Russia LLC Truda 174, 454080 Chelyabinsk. Tel.: +7 (351) 214 5566 / 8811 / 1516

East Africa


Saikang Medical Solutions Ltd. Viraj Complex Unit 3, Mombasa Road, P.O. Box 14262-00800, Nairobi. Kenya. Tel.: +254 733 704 690


Middle East


Saikang Medical Supplies LLC X22-S03 England Cluster, International City, Dubai. United Arab Emirates. Tel.: +20 128 8962497

Healthward uk



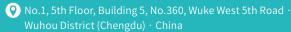

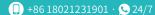

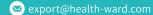

(www.health-ward.com









Chengdu Saikang Medical Equipment Co., Ltd. STOCK CODE: 870098

Company video

General catalogue