

SKB041-10

PATIENT TRANSPORTATION TROLLEY

Registration code: SKM-F



Healthward Tuk

Healthier your life

SKB041-10

FEATURES AND HIGHLIGHTS

Technical parameters

External size (LxW)	1880x620 mm
Table height adjustment	560-890 mm
Back plate adjustment	0-75°
Knee-rest adjustment	0-40°
Revrese trendelenburg	0-18°
Trendelenburg	0-18°

Technical configuration

✓ Crank system	1 set
☑ Back plate air spring	1 set
☑ Plastic side guardrail	1 set
☑ 8" castors with central locking system	4 pcs
☑ Directional wheel	1 pc
☑ Directional pedal	2 pcs
☑ IV pole	1 set
✓ IV pole prevision	4 pcs
Oxygen bottle holder	1 pc
✓ Transfer mattress	1 pc




















Guardrail

Side-rail can be Fixed at Horizontal Level Fixable at horizontal level, where patient's arm can be placed on for the administration of IV or other treatment.

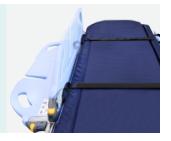

One direction moving

Proper procedure needs to be exercised to safely transfer the patient from the stretcher to the bed. An innovative side-rail design allows medical personnel to carry out this action efficiently and with greater ease. Side-rail can be laid out on a bed positioned next to the stretcher to act as a transfer board, creating zero gap for easy and quick transfer.



Degree indicater

Scale on guardrail, shows the backrest degree.



T-shaped gap

T-shaped design allows tubes to pass through side-rail for improved fluid drainage.



Hydraulic system

American brand hydraulic pump, providing strong support for the functions of height adjustment, trendelenburg/reversed trendelenburg.



Bumper

Equipped with bumper at four



Foldable handle

Equipped with foldable handle.



Center Wheel

Can be engaged to facilitate direction control, enabling users to navigate the stretcher easily in hectic and dynamic environments.



Caster

8" centrolled castors.



I.V pole

Full SS IV pole, bearing 15kg.



Mattress

Wator proof mattress.

















Design

Research

Manufacturing

Maintenance

Service



USA

Saikang Medical Technology, Inc. 17800 Castleton St Ste 180, City of Industry, 91748, California. Tel.: +1 949648 2776

United Kingdom

UK Healthward International Ltd.Dept. 111, 196 High Road Wood
Green, N22 8HH, London.
Tel.: +44 2032399738

Russia

Saikang Medical Russia LLC Truda 174, 454080 Chelyabinsk. Tel.: +7 (351) 2145566 / 8811 / 1516

East Africa

Saikang Medical Solutions Ltd. Viraj Complex Unit 3, Mombasa Road, P.O. Box 14262-00800, Nairobi. Kenya. Tel.: +254733704690

Middle East

Saikang Medical Supplies LLC X22-S03 England Cluster, International City, Dubai. United Arab Emirates. Tel.: +20 128 8962497

Healthward uk

Chengdu Saikang Medical Equipment Co., Ltd. STOCK CODE: 870098

- No.1, 5th Floor, Building 5, No.360, Wuke West 5th Road
- +86 18021231901 · \$\sum_2\$ 24/7
- export@health-ward.com
- www.health-ward.com







Company video

General catalogue