

SK-LLY52C

LED OPERATION LAMP





SK-LLY52C

FEATURES AND HIGHLIGHTS

The LED surgical shadowless lamp is composed of multiple lamp bulb. The bulb has a longer life than the traditional shadowless lamp. It can be round or petal-shaped. It is fixed on the balance arm suspension system. What's more, it is stable in positioning and can be moved vertically or cyclically to meet different heights during surgery.



Technical parameters

Diameter of lighthead (mm):	520
Illumination at 1m (Lux):	160000
Brightness control range:	0-100% stepless adjustment
Color temperature (Kelvin):	3000-6700
Color rendering index (CRI):	≥98
 Temperature rise (Surgeon head)(°C): 	≤0.5
 Temperature rise (Surgical field)(°C): 	≤1.5
Total irradiance Ee (W/m2):	≤500
 Radiant power Ee/Ec (mW/m2lux): 	≪4
Effective focusing depth (cm):	130
Lighting spot size (cm):	10-35
• L1+L2 depth of illumination (cm):	≥100
Bulb type:	LED

Bulb brand:	OSRAM
Bulb quantity:	54
■ Bulb power (w):	3
Average bulb life (hours):	≥60000

Technical configuration

✓ Rotating arm	1 set
☑ Balance arm	1 set
☑ C-arm	1 set
☑ Handle	1 pcs



Material technology

The lampshade shell adopts a thin, light and streamlined design, and ABS is integrally formed to achieve the best laminar flow effect and meet the requirements of modern clean laminar flow operating rooms.

The surface of the lamp arm adopts the high-voltage electrostatic spraying process, and the imported environmentally friendly antibacterial plastic powder that used to ensures that the product meets the requirements of surgical hygiene.



Bulb


Using the original imported LED lamp beads from Osram, Germany, the average service life is as high as 60,000 hours, which is 40 times longer life than halogen lamps.

Under the same brightness, the energy consumption of LED is only 1/10 of that of ordinary incandescent lamps and 1/2 of that of halogen lamps.



Color temperature&illumination

Using a new type of LED cold light source, the illuminance can reach 160,000 Lux, which can realize stepless dimming and non-multi-speed adjustment; the color temperature is in the range of 3000K-6700K, and the stepless dimming is non-multi-step adjustment. While meeting the high illuminance, adjust the lighting parameters according to the needs of different surgeons, so that the surgeon's perception of light is soft and no glint.

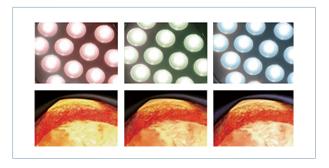

Color rendering index


The color rendering index reaches 98, which truly reflects the color of human tissues. It is suitable for various surgical scenes and greatly reduces the visual fatigue of medical staff due to long-term operations.






Temperature rise (surgical field)

LED cold light source has no infrared radiation, coupled with our nano-coated radiator, created our excellent heat dissipation effect. It solves the problem of the doctor's head temperature being too high during the operation, and creates a more comfortable operation environment for the doctor; it also effectively avoids the loss of water and infection of the surgical wound, which is conducive to the rapid recovery of the patient's wound after the operation.

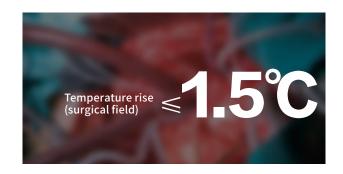

Reflective bowl design

The multi-focus design of the reflective bowl is adopted, and 3 luminous bodies are hidden in each reflective bowl. The light source in the reflective bowl is just rightly converged to the surgical area, forming a good shadowless effect.

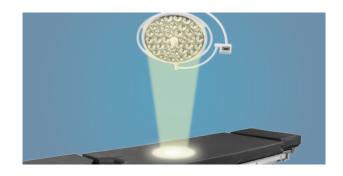

Operating handle

The detachable disinfection handle is made of high-strength ABS, which can be disinfected and sterilized at a high temperature above 135°, ensuring the durability and reliability of the handle.



Other functions

With one-key cavity mirror mode, can easily switch between standard mode and cavity mirror mode.













Research



Manufacturing



Maintenance





Saikang Medical Technology, Inc. 17800 Castleton St Ste 180, City of Industry, 91748, California. Tel.: +1 949648 2776

United Kingdom

 ${\tt UK\, Healthward\, International\, Ltd.}$ Dept. 111, 196 High Road Wood Green, N22 8HH, London. Tel.: +44 2032399738

Russia

Saikang Medical Russia LLC Truda 174, 454080 Chelyabinsk. Tel.: +7 (351) 214 5566 / 8811 / 1516

East Africa

Saikang Medical Solutions Ltd. Viraj Complex Unit 3, Mombasa Road, P.O. Box 14262-00800, Nairobi. Kenya. Tel.: +254 733 704 690

Middle East

Saikang Medical Supplies LLC X22-S03 England Cluster, International City, Dubai. United Arab Emirates Tel.: +20 128 8962497



Jiangsu Saikang Medical Equipment Co., Ltd. STOCK CODE: 870098

No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China

export@sksurgitech.com

www.sksurgitech.com















General catalogue